

5-part: laser + impedance + colourimetric haematology

Three discrete measurement channels combine to provide a true 5-part differential, with White Blood Cell, Red Blood Cell and Platelet parameter results from just a 15µL sample. In only 60 seconds get results you can trust, with Laser flow cytometry, impedance counts and colourimetric detection technologies to detect not only white blood cell volume, but complexity and granularity obtained by a triple angle laser scatter. While accurate and fast, the Element HT5 remains very simple to use with a large colour touch screen, intuitive navigation, simplified data entry and easy reagent management.





www.heska.com.au

Combination of Technologies

	Counting Channels		Technology	Parameters
	1	Diff Channel	Flow Cytometry	LYM, MON, NEU, EOS
	2	RBC/PLT Channel	Impedance	RBC/PLT
	3	Baso/HGB Channel	Colourimetry	BAS, WBC, HGB

Features:

Consolidated results and reporting

• Combine results using the HeskaView Integrated Software® program or transmit to your PMS. Rapidly retrieve patient worklist with bidirectional communication.

Multi-species availability

• Dog, Cat, Horse, Cow, Ferret, Goat, Llama, Monkey, Mouse, Giant Panda, Red Panda, Pig, Rabbit, Rat and Sheep

Test Parameters

White Blood Cell	Red Blood Cell	
Total White Blood Cell (WBC)	Total Red Blood Cell (RBC)	
Neutrophil Number (NEU#)	Haemoglobin (HGB)	
Neutrophil Percent (NEU%)	Haematocrit (HCT)	
Lymphocyte Number (LYM#)	Mean Cell Volume (MCV)	
Lymphocyte Percent (LYM%)	Mean Cell Haemoglobin (MCH)	
Monocyte Number (MON#)	Mean Cell Haemoglobin Conc. (MCHC)	
Monocyte Percent (MON%)	Red Blood Cell Distribution (RDW-CV)	
Eosinophil Number (EOS#)	—	
Eosinophil Percent (EOS%)	Platelet	
Basophil Number (BAS#)	Total Platelet (PLT)	
Basophil Percent (BAS%)	Mean Platelet Volume (MPV)	

Results and

Easy to use

- Web-style links for simple navigation
- Reporting
- Simple data management Rapidly retrieve trending data from patient record

Efficient

- Assign, unassign or delete multiple items at once
- Species-specific reference ranges and colour-coded flags make it easy to quickly evaluate patient data

www.heska.com.au 1300 437 522

©2019 Heska Corporation. All Rights Reserved. HESKA, Element HT5 and HeskaView Integrated Software are registered trademarks of Heska Corporation. AU190208-001

Small Sample Size of 15µL




Easy reagent management

• Manage reagents via bar-code scanner and on-board menu functions. Internal reagent storage saves space and reduces clutter


Intuitive navigation


• Simplifed data entry and patient test review with numeric results, histograms and scatter plots



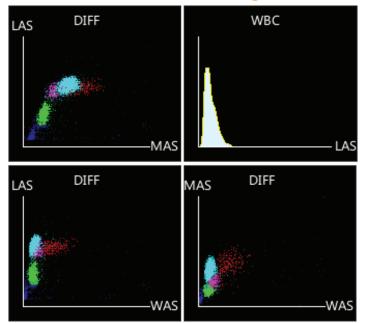

RBC /PLT Histograms

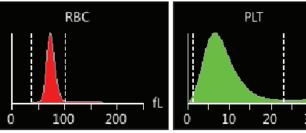

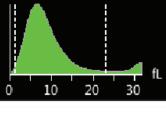



WBC Scatter Plots & Histograms